

# Attends<sup>®</sup> Slip Regular

## Question:

When does the bladder weakness occur?

What happens?

**More comfort. More dryness. More security.**

The new, improved Attends Slip Regular has a powerful Quick-Dry absorbent zone to provide all-round dryness, helping you to feel even more comfortable. The active core has Attends tried-and-tested Quick-Dry system as well as a third absorbent core. Leakage barriers ensure maximum side protection.

Unpleasant odours? Not any more. The Attends Odour Protection™ system binds the ammonia found in urine and locks it away inside the slip.

Breathable non-woven wings allow air to circulate. Leg elastics ensure the perfect fit. With its 4 fixation tapes, the Attends Slip Regular is easy to fit and comfortable to wear.

## Answer:

The client has no control over bladder and/or bowel function.


Heavy losses of urine and/or faeces, at any time during the day or night.



		9 X-Small	9 Small	8 Medium	8 Large	10 Extra Large
Size / waist size	cm	40-70	50-80	80-110	100-150	150-175
Packcount	Pcs	1 x 28	1 x 28	1 x 28	1 x 28	1 x 14
Cascount	Pcks	6 x 28	6 x 28	4 x 28	4 x 28	4 x 14
Cases per pallet	Case	16 x 6 x 28	16 x 6 x 28	24 x 4 x 28	16 x 4 x 28	28 x 4 x 14
Length	cm	58	71	79	96	103
Weight	g	78	86	85	104	179
Absorption capacity (ISO)	ml	1500	1700	1,700	2,100	3600
Quick-Dry	Material	Super-absorbent fluff	Super-absorbent fluff	Super-absorbent fluff	Super-absorbent fluff	Super-absorbent fluff
Top sheet	Material	Non-woven	Non-woven	Non-woven	Non-woven	Non-woven
Back sheet with wetness indicator		Polyethylene	Polyethylene	Polyethylene	Polyethylene	Polyethylene
Leg Elastics/Leakage Barrier		Synthetic	Synthetic	Synthetic	Synthetic	Synthetic, Lycra
Fixation tapes		Yes	Yes	Yes	Yes	Yes
Quick-Dry technology		Yes	Yes	Yes	Yes	-
Odour Protection™		Yes	Yes	Yes	Yes	Yes

# Attends<sup>®</sup> Slip Regular


## Quick-Dry<sup>®</sup> absorbent zone

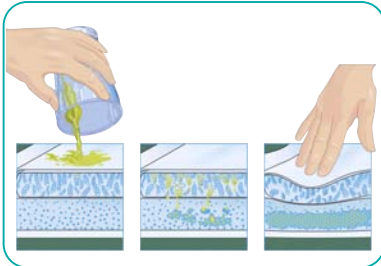

Thanks to its Quick-Dry<sup>®</sup> absorbent zone\*, the Attends Slip ensures allround, super-dry comfort. In addition to the tried-and-tested Attends Quick-Dry<sup>®</sup> system, this innovative technology also incorporates a highly-concentrated third core with super-absorber. Attends Quick-Dry<sup>®</sup> is cutting-edge technology consisting of an acquisition patch in the centre of the slip and its underlying core.

First, Attends Quick-Dry<sup>®</sup> allows the fluids to spread out across the non-woven layer, thus increasing absorption considerably. Then, curly fibres below the nonwoven layer open up to absorb the fluids quickly and lead them to the special internal core, where they are turned into a gel. Afterwards, the fibres close up to ensure that no fluids can leak back to the surface, thus enabling Attends Quick-Dry<sup>®</sup> to reduce the risk of leakage.

## Odour Protection<sup>™</sup>

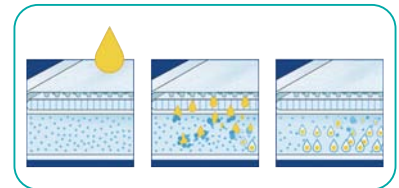

Attends Quick-Dry<sup>®</sup> can also absorb several gushes and lock them away quickly and securely, ensuring optimum dryness for your skin.

The main cause of unpleasant odours is the ammonia in urine, which builds up due to normal bacteria found in the body. The Attends core binds the ammonia found in urine and locks it away, reducing the risk of unpleasant smells.

Binding of ammonia results in a lower pH, more towards the natural pH of the skin. No perfumes or other additives are used.